

Gusta II NC PTCA Balloon Dilatation Catheter

Non-compliant Balloon Dilatation Catheter

Specifications

Balloon Diameter (mm)	Balloon Length(mm)							
	6	8	10	12	15	20	25	30
2.25	DHC2206 □	DHC2208 □	DHC2210 □	DHC2212 □	DHC2215 □	DHC2220 □	DHC2225 □	DHC2230 □
2.5	DHC2506 □	DHC2508 □	DHC2510 □	DHC2512 □	DHC2515 □	DHC2520 □	DHC2525 □	DHC2530 □
2.75	DHC2706 □	DHC2708 □	DHC2710 □	DHC2712 □	DHC2715 □	DHC2720 □	DHC2725 □	DHC2730 □
3.0	DHC3006 □	DHC3008 □	DHC3010 □	DHC3012 □	DHC3015 □	DHC3020 □	DHC3025 □	DHC3030 □
3.25	DHC3206 □	DHC3208 □	DHC3210 □	DHC3212 □	DHC3215 □	DHC3220 □	DHC3225 □	DHC3230 □
3.5	DHC3506 □	DHC3508 □	DHC3510 □	DHC3512 □	DHC3515 □	DHC3520 □	DHC3525 □	DHC3530 □
3.75	DHC3706 □	DHC3708 □	DHC3710 □	DHC3712 □	DHC3715 □	DHC3720 □	DHC3725 □	DHC3730 □
4.0	DHC4006 □	DHC4008 □	DHC4010 □	DHC4012 □	DHC4015 □	DHC4020 □	DHC4025 □	DHC4030 □
4.25	DHC4206 □	DHC4208 □	DHC4210 □	DHC4212 □	DHC4215 □	DHC4220 □	DHC4225 □	DHC4230 □
4.5	DHC4506 □	DHC4508 □	DHC4510 □	DHC4512 □	DHC4515 □	DHC4520 □	DHC4525 □	DHC4530 □

Compliance Chart

Pressure (atm)	Balloon Diameter(mm)									
	2.25	2.5	2.75	3.0	3.25	3.5	3.75	4.0	4.25	4.5
6	2.16	2.35	2.63	2.81	3.08	3.23	3.56	3.62	4.04	4.12
8	2.19	2.42	2.67	2.88	3.14	3.37	3.63	3.80	4.11	4.27
10	2.22	2.46	2.71	2.94	3.20	3.43	3.69	3.92	4.18	4.40
12(NP)	2.25	2.50	2.75	3.00	3.25	3.50	3.75	4.00	4.25	4.50
14	2.28	2.56	2.79	3.07	3.30	3.58	3.81	4.06	4.32	4.56
16	2.31	2.62	2.83	3.13	3.35	3.64	3.87	4.12	4.38	4.62
18	2.34	2.67	2.88	3.18	3.40	3.70	3.93	4.18	4.44	4.68
20(RBP)*	2.37	2.72	2.91	3.23	3.44	3.77	3.98	4.25	4.50	4.74
22	2.40	2.77	2.94	3.28	3.48	3.84	4.03	4.31	4.56	4.80

\*Rated Burst Pressure.Do not Exceed.

Gusta II PTCA Balloon Dilatation Catheter

Semi-compliant Balloon Dilatation Catheter

Specifications

Balloon Diameter (mm)	Balloon Length(mm)							
	6	8	10	12	15	20	25	30
1.25	DPC1206 □	DPC1208 □	DPC1210 □	DPC1212 □	DPC1215 □	DPC1220 □	-	-
1.5	DPC1506 □	DPC1508 □	DPC1510 □	DPC1512 □	DPC1515 □	DPC1520 □	DPC1525 □	DPC1530 □
1.75	DPC1706 □	DPC1708 □	DPC1710 □	DPC1712 □	DPC1715 □	DPC1720 □	DPC1725 □	DPC1730 □
2.0	DPC2006 □	DPC2008 □	DPC2010 □	DPC2012 □	DPC2015 □	DPC2020 □	DPC2025 □	DPC2030 □
2.25	DPC2206 □	DPC2208 □	DPC2210 □	DPC2212 □	DPC2215 □	DPC2220 □	DPC2225 □	DPC2230 □
2.5	DPC2506 □	DPC2508 □	DPC2510 □	DPC2512 □	DPC2515 □	DPC2520 □	DPC2525 □	DPC2530 □
2.75	DPC2706 □	DPC2708 □	DPC2710 □	DPC2712 □	DPC2715 □	DPC2720 □	DPC2725 □	DPC2730 □
3.0	DPC3006 □	DPC3008 □	DPC3010 □	DPC3012 □	DPC3015 □	DPC3020 □	DPC3025 □	DPC3030 □
3.25	DPC3206 □	-	DPC3210 □	DPC3212 □	DPC3215 □	DPC3220 □	DPC3225 □	DPC3230 □
3.5	DPC3506 □	-	DPC3510 □	DPC3512 □	DPC3515 □	DPC3520 □	DPC3525 □	DPC3530 □
3.75	DPC3706 □	-	DPC3710 □	DPC3712 □	DPC3715 □	DPC3720 □	DPC3725 □	DPC3730 □
4.0	DPC4006 □	-	DPC4010 □	DPC4012 □	DPC4015 □	DPC4020 □	DPC4025 □	DPC4030 □
4.5	-	-	-	-	DPC4515 □	DPC4520 □	-	-
5.0	-	-	-	-	DPC5015 □	DPC5020 □	-	-

Compliance Chart

Pressure (atm)	Balloon Diameter(mm)														
	1.25	1.5	1.75	2.0	2.25	2.5	2.75	3.0	3.25	3.5	3.75	4.0	4.5	5.0	
4	1.17	1.43	1.68	1.91	2.15	2.39	2.60	2.87	3.09	3.33	3.68	3.81	4.31	4.79	
6(NP)	1.25	1.50	1.75	2.00	2.25	2.50	2.75	3.00	3.25	3.50	3.75	4.00	4.50	5.00	
8	1.31	1.57	1.80	2.07	2.32	2.58	2.81	3.12	3.36	3.62	3.83	4.09	4.61	5.11	
10	1.36	1.63	1.86	2.12	2.37	2.63	2.87	3.20	3.45	3.72	3.91	4.20	4.71	5.21	
12	1.39	1.67	1.92	2.16	2.41	2.69	2.93	3.27	3.52	3.80	4.01	4.29	4.82	5.32	
14	1.43*	1.71*	1.97	2.20	2.45	2.75	2.97	3.34	3.58	3.89	4.12	4.39	4.93	5.43	
16(RBP*)	1.50	1.75	2.03	2.24	2.50	2.81	3.05	3.41	3.63	3.98	4.23	4.51	5.04	5.56	
18	1.53	1.78	2.11	2.29	2.55	2.88	3.10	3.50	3.68	4.06	4.39	4.61	5.15	5.72	

\*Rated Burst Pressure.Do not Exceed.The part marked in blue is the RBP of balloon.



Demax Medical Australia Pty Ltd.  
B17/148 Old Pittwater Rd Brookvale 2100

www.demax.group  
info@demaxmedical.com



DIO MEDICAL  
76, Sangamsan-ro, Mapo-gu, Seoul, Korea

www.diomedi.co.kr  
info@diomedi.co.kr



Gusta


PTCA Balloon Dilatation Catheter

Semi-compliant/Non-compliant





# ARTERIES

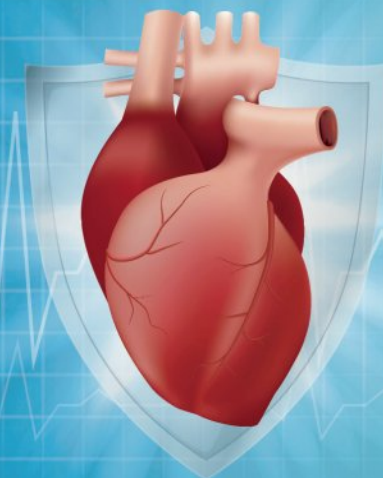

## Gusta PTCA Balloon Dilatation Catheter

For the dilatation of coronary stenosis or intravascular stent




# Gusta PTCA



- ① **Tapered tip design**  
Small lesion entry profile (0.016" = 0.4mm)  
Tightly wrapped to the guidewire, avoiding fish-mouth phenomenon while ensuring traceability and crossability
- ② **Embedded platinum-iridium alloy marker**  
Reduces outer diameter loss, improves crossability and ensures precise positioning
- ③ **Selective hydrophilic coating\***  
Provides lasting lubrication for tracking and lesion-crossing before balloon dilatation  
The uncoated portion reduces balloon sliding after balloon dilatation (\*Hydrophilic coating patented)
- ④ **Smooth guidewire exit port with a smaller outer diameter**  
For easier delivery and reduced intravascular resistance
- ⑤ **Metal support shaft (hypotube) that tapers toward the distal end**  
Balancing support and flexibility to ensure the effective transfer of pushing force without bending

## Gusta PTCA Balloon Dilatation Catheter

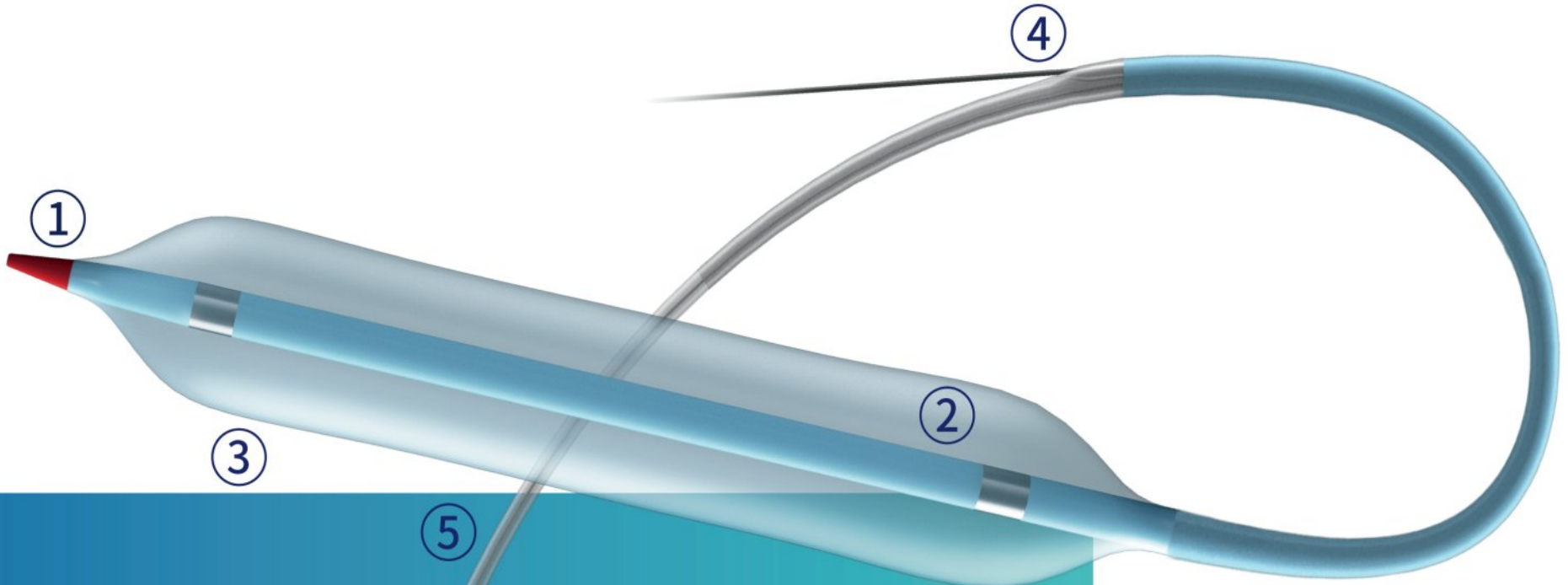

Naming System (DPC2515)		
DPC	25	15
Specification	Diameter	Length


POT: Gusta NC  
4.0\*12mm

Pre-dilatation: Gusta  
1.25\*15mm


Kissing: Gusta NC  
3.5\*12mm, 3.0\*12mm

RCA case: Gusta  
2.5\*20mm

### Application cases

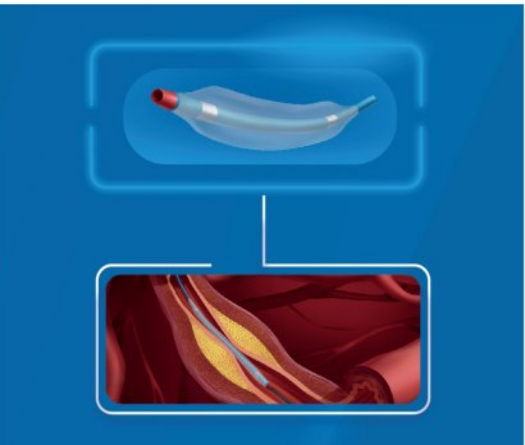

- ⊕ Excellent crossability and good pushability
- ⊕ Accurate compliance and safe burst pressure ensure ideal dilatation outcome
- ⊕ A broad range of specifications, allows for refined manipulation when treating complex lesions

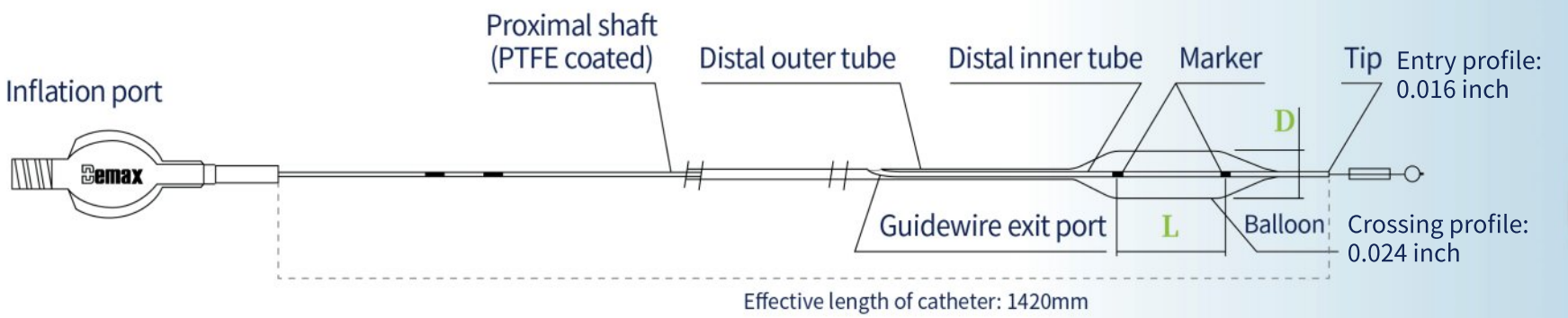

Note: Balloon dilatation catheters with 1.25 and 1.5 mm balloon diameters, only one marker.